




Medvance[®]

Silicone Foam

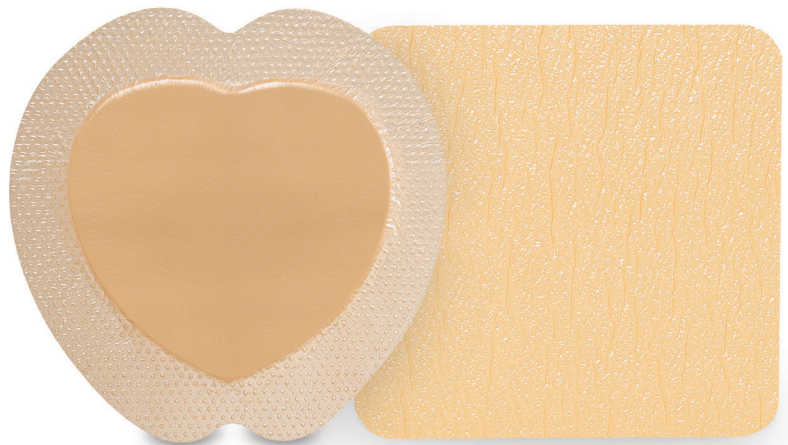

Improving patient care and enhancing the quality of people's lives



MedVANCE[®] Silicone Foam

Product Description:

MedVance[®] Silicone Foam Dressing is an absorbent self-adhering dressing consisting of a polyurethane (PU) foam pad, a soft, hypoallergenic silicone adhesive wound contact layer, and a high MVTR transparent film top layer.

The multilayer structure can absorb large amounts of wound exudate, and the soft silicone adhesive acts as a gentle wound contact layer that seals the wound edges but does not stick to the wound bed, minimizing pain and medical adhesive related skin damage (MARS!) when changing dressings.

Breathable PU film allows the passage of oxygen and water vapor, extending wear time and ensuring a moist wound healing environment.

MedVance[®] Silicone Foam Dressings are available with and without an adhesive border, and as a conformable exudate transfer dressing.

Features and Benefits:

- Soft silicone adhesive adheres to the skin around the wound but not the wound.
- Easy release adhesive minimizes pain and adhesive associated skin damage on removal.
- Soft foam flexes and conforms to body contours for maximum patient comfort.
- High-capacity foam extends wear time and combats maceration.
- Low-friction transparent film reduces shear and protects from contamination and moisture-associated skin damage while maintaining a moist wound healing environment.
- Unique highly extensible wing format remains secure over mobile locations.




Indications for use:

MedVance[®] Silicone Foam Dressing is indicated for the management of moderate to highly exuding full or partial thickness wounds such as:

- Pressure injuries
- Venous stasis ulcers
- Diabetic foot ulcers
- Graft donor sites, abrasions, and other wounds healing by secondary intention

5 Layers of protection

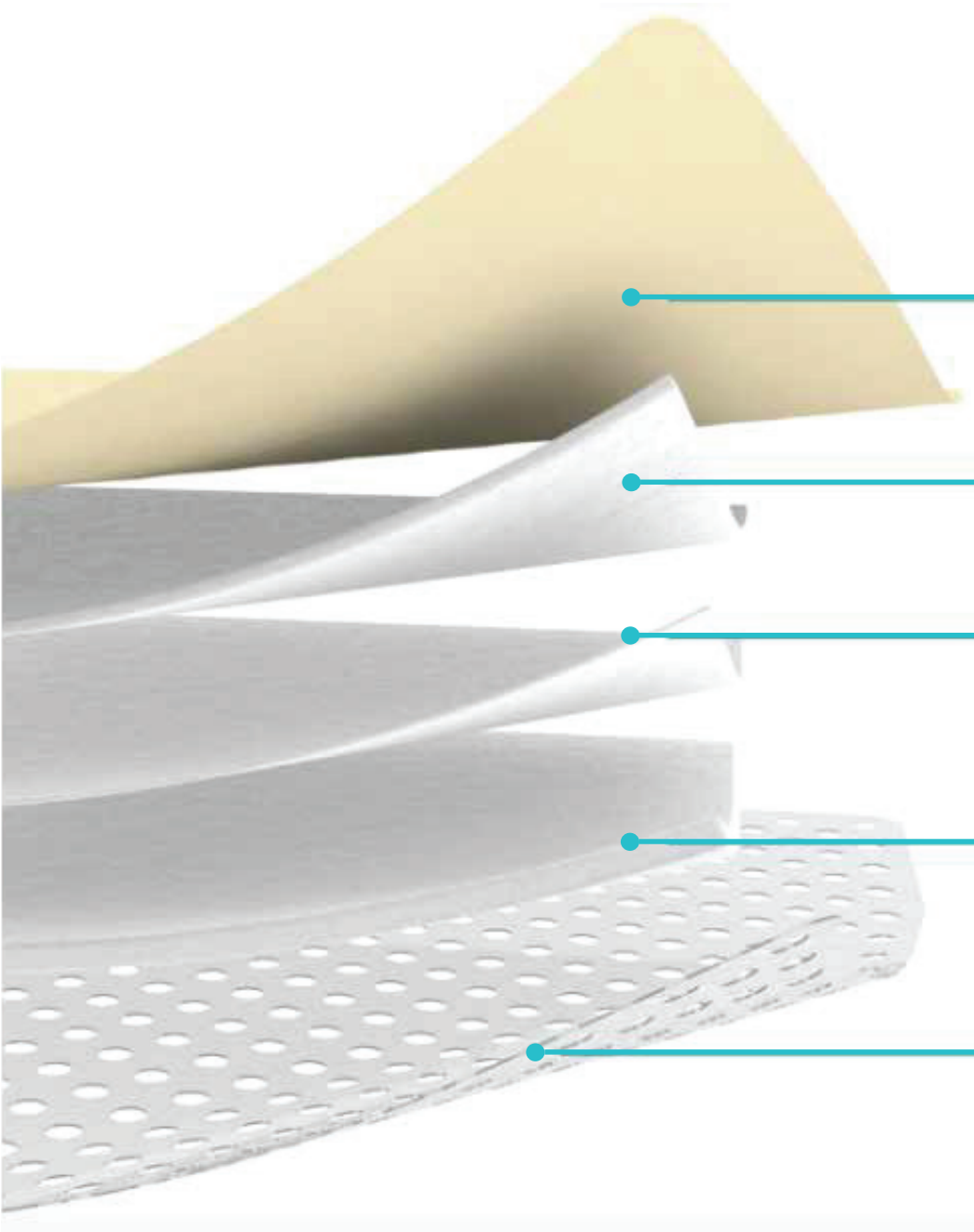

Protective Backing

- Fluid and microbial barrier
- Conformable with low coefficient of friction
- High moisture vapor transmission rate

Super Absorbent Core

- High fluid absorption capacity
- Retains fluid under compression
- Extends dressing wear time

Nonwoven Distribution Layer


- Distributes fluid over entire dressing
- Maximizes use of the absorbent core
- Extends wear time

Polyurethane Foam

- Broad spectrum
- Retains potency for a minimum of 7 days
- Soft foam comfort and pressure relief

Silicone Wound Contact Layer

- Gentle and secure hypoallergenic silicone adhesive
- Adheres to skin but not to fragile wound tissue
- Provides pain-free, atraumatic removal

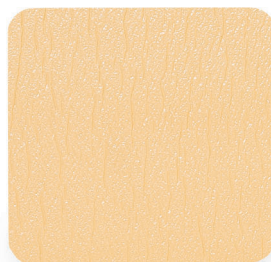

Product Availability

MedVance Silicone Foam Dressing

Product Code	Size	Packaging
V3000202	2" x 2"	5 pcs/box 80 boxes/case
V3000303	3" x 3"	5 pcs/box 80 boxes/case
V3000404	4" x 4"	5 pcs/box 80 boxes/case
V3000606	6" x 6"	5 pcs/box 80 boxes/case
V3000405	4" x 5"	5 pcs/box 80 boxes/case
V3000408	4" x 8"	5 pcs/box 80 boxes/case
V3000608	6" x 8"	5 pcs/box 40 boxes/case
V3000808	8" x 8"	5 pcs/box 40 boxes/case

MedVance Silicone Foam Dressing With Border

Product Code	Size	Packaging
V3010202	2" x 2"	5 pcs/box 80 boxes/case
V3010205	2" x 5"	5 pcs/box 80 boxes/case
V3010303	3" x 3"	5 pcs/box 80 boxes/case
V3010404	4" x 4"	5 pcs/box 80 boxes/case
V3010408	4" x 8"	5 pcs/box 80 boxes/case
V3010412	4" x 12"	5 pcs/box 80 boxes/case
V3010606	6" x 6"	5 pcs/box 80 boxes/case
V3010608	6" x 8"	5 pcs/box 80 boxes/case
V3010707	7" x 7"	5 pcs/box 80 boxes/case
V3210508 Heel	5" x 8"	5 pcs/box 80 boxes/case
V3110707 Sacral sm.	7" x 7"	5 pcs/box 80 boxes/case
V3110909 Sacral lr.	9" x 9"	5 pcs/box 40 boxes/case

MedVance Silicone Foam Transfer Dressing

Product Code	Size	Packaging
V3020820	8" x 20"	5 pcs/box 20 boxes/case

MedWay Inc.

1650 Horizon Parkway, Suite 450
 Buford, GA 30518
 (888) 5MEDWAY (633929)
 sales@medwayus.com
 www.medwayus.com

mv.br.3.1-v1