



MedAir[®]

Elite




Improving patient care and enhancing the quality of people's lives



The control and automation you need in an alternating pressure mattress.

1. To-the-pound firmness adjustability
2. Timed alternating air cells to change pressure points at different frequencies.
3. Panel Lock to avoid accidental adjustments

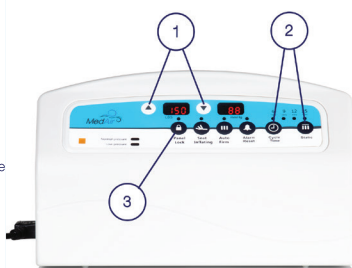

Product Description:

The MedAir Elite Dual Therapy Group II mattress system delivers immersive alternating pressure and low air loss micro-climate management in a versatile and cost-effective design.

The fully customizable digital pump adjusts air flow and pressure redistribution to individual patient requirements while diagnostic alarm settings inform caregivers when patient attention is required.

Auto-return function automatically shifts from static to alternating-pressure mode to ensure maximum pressure redistribution and relief.

The 8" deep cell-on-cell configuration and pillow function provides immersive pressure relief and prevents bottoming out while providing 24-hour power outage protection.

Adjustable low air loss settings deliver low interface pressures and regulate temperature and moisture microclimate to reduce the most common factors leading to pressure injuries and the impediments to healing optimizing patient comfort, wound healing and pressure injury prevention.

Indications for Use:

The MedAir Elite Dual Therapy combination pressure mattress is intended for the management or prevention of pressure injuries as well as the reduction of moisture and temperature at the mattress/skin interface.

Features and Benefits:

- Digital controls
- Intuitive 20" digital LED screen displays weight and pressure
- Automatic pressure control maximizes offloading
- Ultra-silent design delivers noise levels less than 20 dB
- Double pump cores provide maximum air output
- Audio and visual alarms
- Aluminum pump core case accelerates heat dissipation and minimizes vibration
- Static-mode alarm ensures return to alternating pressure mode
- Multiple settings suit every need, including lock mode, sitting mode, remote control, auto firm and more
- Remote control
- Waterproof quilted top cover provides four-way stretch and zipper closure

Product Details:

Product code: MA-M4P8-US

HCPCS: E0277

Specifications:

- Dimensions: 80"(L) x 36"(W) x 8"(H)
- Pump airflow: 8 LPM
- Cycle times available: 6, 9, 12 or 25 minutes
- Weight capacity: 450 lbs

MedWay Inc.

1650 Horizon Parkway, Suite 450
 Buford, GA 30518
 (888) 5MEDWAY (633929)
 sales@medwayus.com
 www.medwayus.com

MA-M4P8-US